

## Application News


MALDI-2201

High Mass MALDI-TOF Mass Spectrometry

### Characterization of intact monoclonal Antibody/Antigen complexes by the OmegaToF MS

#### ■ Introduction

The ability of MALDI-TOF mass spectrometry to accurately perform protein detection is well recognized in the life science field. The OmegaToF MALDI mass spectrometer has made it now possible to analyze intact proteins and protein interactions up to 1500 kDa. This extended mass range detection capability allows for the characterization of most non-covalent protein complexes, including the antigen/antibody interactions essential for the development of protein therapeutics.

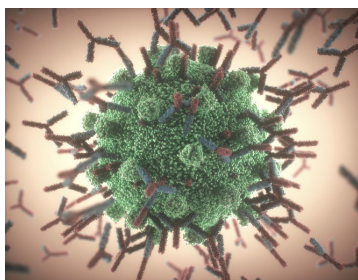

The extended mass range of the OmegaToF gives direct access to the characterization of complexes formed between therapeutic monoclonal antibodies (mAbs) and their antigen targets. The analysis becomes routine as the samples are analyzed directly in their original buffers, with the sensitivity and robustness of standard MALDI ToF mass spectrometric analysis.

Here, we demonstrate the capability of the OmegaToF mass spectrometer (OmegaTOF, Fig. 1) to perform high-mass detection of intact antigen/mAb protein complexes.



**Figure 1:** OmegaToF MALDI Mass Spectrometer

### ■ Sample Preparation and Measurement Conditions

Antigen and antibody samples were purchased from Sino Biological. These were prepared at a final concentration of 1 pmol/μL each in the final mixture. Fifty shots were accumulated per profile (200 profiles per spectrum). The mass spectra were recorded using the average masses.

For cross-link experiments, 9 μL of the antibody/antigen sample prepared previously is mixed with 1 μL of R200 cross-linking reagent [CovalX R200] to reach a final concentration of 1 mg/mL. The sample is incubated 180 minutes before direct mass spectrometric analysis in the OmegaToF.

### ■ Antibody : Antigen Complex

To demonstrate the OmegaTOF sensitivity, the antigen/antibody complex spectra are presented. Fig. 3 shows the mass spectrum of the antigen/antibody mixture. The singly-charged ion of the antigen (approx. 107kDa) and the singly-charged ion of the antibody (approx. 143kDa) were observed with good signal-to-noise ratio.

Fig. 4 shows the mass spectrum of the same antigen/antibody mixture analyzed in Fig. 3 after cross-linking. The singly-charged ion of the antigen/Antibody complex (approx. 258kDa, stoichiometry 1:1) and the singly-charged ion of the antibody (approx. 148kDa) were observed. The singly-charged ion of the antigen observed in Fig. 3 is not observed after crosslinking due to the antibody binding.

Fig. 5 shows an overlay of the control (Fig. 3) and cross-link spectra (Fig. 4). The singly-charged ion (approx. 248 kDa) of the antigen/antibody complex with stoichiometry 1:1 is clearly observed as the primary signal after crosslinking. The stoichiometry of 2:1 is not observed for this experiment.

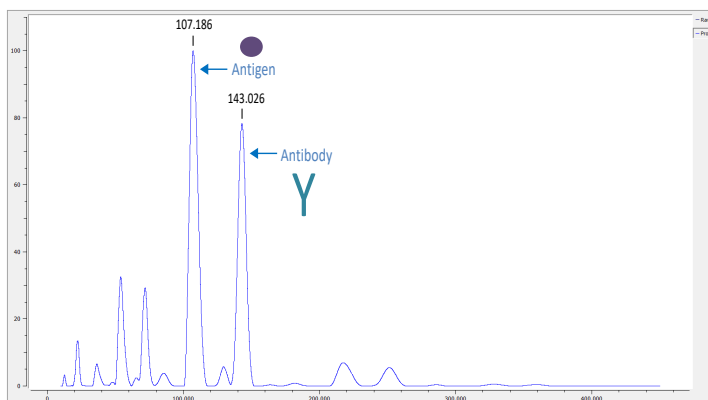


Figure 3: OmegaToF spectrum of antibody/antigen mixture.

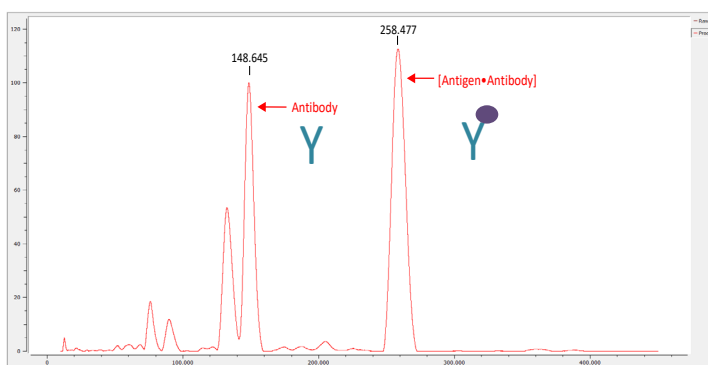


Figure 4: Detection of antibody-antigen complex after cross-linking.

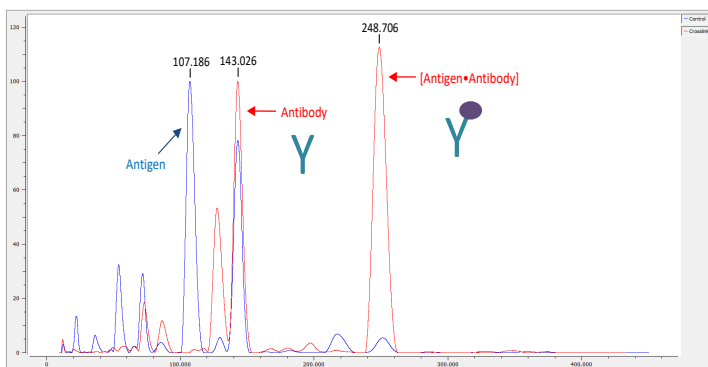



Figure 5: Overlay (control, blue; cross-link, red) of antibody-antigen complex.

## ■ Two Antibodies Binding One Antigen

This technology can also be used to characterize whether two different monoclonal antibodies are able to bind onto the antigen protein at the same time. If they are able to bind at the same time, then this signifies that they are binding onto different epitope regions of the antigen protein.

Fig. 6 shows an overlay of two mass spectra containing two different monoclonal antibodies against the same target antigen protein. In blue is the control spectra of the two monoclonal antibodies and the antigen mixed together and analyzed together. In the blue (control) spectrum, the singly-charged ions of Antibody 1 or Antibody 2 (approx. 147 kDa) and the singly charged ion of the Antigen (approx. 87 kDa) were observed.

The red spectrum is the signal of the crosslinked material from the same control (blue) sample. The detection of a complex binding both monoclonal antibodies (approx 384kDa) and the antigen together at the same time provides insight that the two monoclonal antibodies are binding different epitope regions on the antigen protein.

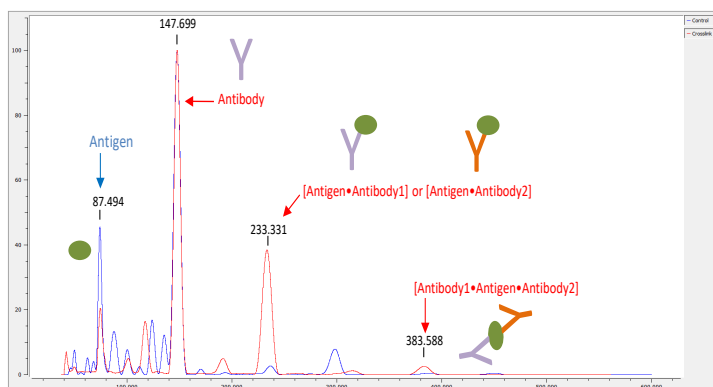

**Figure 6:** Overlay (control, blue; crosslink, red) of mixture containing two different monoclonal antibodies (Antibody 1 and Antibody 2) and target antigen.



First Edition: February 2022

**SHIMADZU Corporation**  
www.shimadzu.com/an/

**SHIMADZU SCIENTIFIC INSTRUMENTS**  
7102 Riverwood Drive, Columbia, MD 21046, USA  
Phone: 800-477-1227/410-381-1227, Fax: 410-381-1222  
URL: www.ssi.shimadzu.com

For Research Use Only. Not for use in diagnostic procedure.  
This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. Shimadzu disclaims any proprietary interest in trademarks and trade names used in this publication other than its own. See <http://www.shimadzu.com/about/trademarks/index.html> for details.

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.