

Accurate Mass LC-MS/MS Profiling of Synthetic Cannabinoids

IMSC 2012 PTu-089

Ben Figard¹; Christopher L. Pennington²; Jeff Dahl³;
Susan Leonard⁴; Jorge Smith¹,

¹Shimadzu Scientific Instruments, Houston, TX;

²Rice University, Houston, TX;

³Shimadzu Scientific Instruments, Columbia, MD;

⁴Shimadzu Scientific Instruments, Marlborough, MA.

Accurate Mass LC-MS/MS Profiling of Synthetic Cannabinoids

Introduction

New synthetic designer drugs of abuse are being produced as readily available consumer products such as bath salts and incense. The uncontrolled and unregulated synthesis of these drugs causes numerous variations in formula, structure, and stereochemistry. These variations cause the

analysis of the compounds of interest to be very difficult and extremely time consuming. In this study we describe a new simple and efficient method by which these new synthetic drugs of abuse can be rapidly identified and quantified.

Materials and Method

The samples received were herbal in nature. It was decided that aliquots of each sample would be extracted in acidic

aqueous, basic aqueous and various organic solvents. The resulting extracts were then analyzed by LC-MS.

Sample Preparation

For each extraction 224 – 248 mg of sample was weighed into a 20 mL scintillation vial, 5 mL of extraction solvent was added, samples were vortexed and sonicated for 60 minutes. Following sonication the extraction solvent was removed to another vial and (with the exception of aqueous solvents) evaporated to dryness under nitrogen. Solvents used for extraction were as follows: CHCl₃,

Acetone, IPA, MeOH, ACN, 0.1 % NH₄OH in water and 0.1 % formic acid in water. For LC-MS analysis samples were reconstituted with the addition of 2 mL of 0.1% formic acid in water followed by 3 mL of acetonitrile. Following solvent addition samples were vortexed, sonicated and analyzed by LC-MS. Aqueous extracts were analyzed as is.

Chromatography

Instrument : Shimadzu XR system
 Column : Supelco Ascentis Express C-18 150 mm × 2.1 mm, 2.7 μm
 Mobile Phase : A: 0.1% Formic Acid in Water
 B: 0.1% Formic Acid in ACN
 Flow Rate : 0.4 mL/min
 Injection Volume : 25 μL

Gradient Program:

Gradient A		Gradient B	
Time (min)	%B	Time (min)	%B
0	5	0	30
35	20	45	90
49	90	54	90
54	90	55	30
55	5	60	30
60	5		

Accurate Mass LC-MS/MS Profiling of Synthetic Cannabinoids

Mass Spectrometry

Instrument : Shimadzu LCMS-IT-ToF Mass Spectrometer

Ionization : ESI

Polarity : Positive


Mass Range: 100 – 500 m/z

Data Analysis


ESI-MS accurate mass data was used to predict chemical formulas which were searched in databases such as ChemSpider (<http://www.chemspider.com/>) and MassBank (<http://www.massbank.jp>).

ESI-MS accurate mass data and MS/MS data were also searched using METLIN (<http://metlin.scripps.edu/>). MassFrontier 5.0 (ThermoFisher) was used to help evaluate MS/MS data and propose chemical structures.

Results

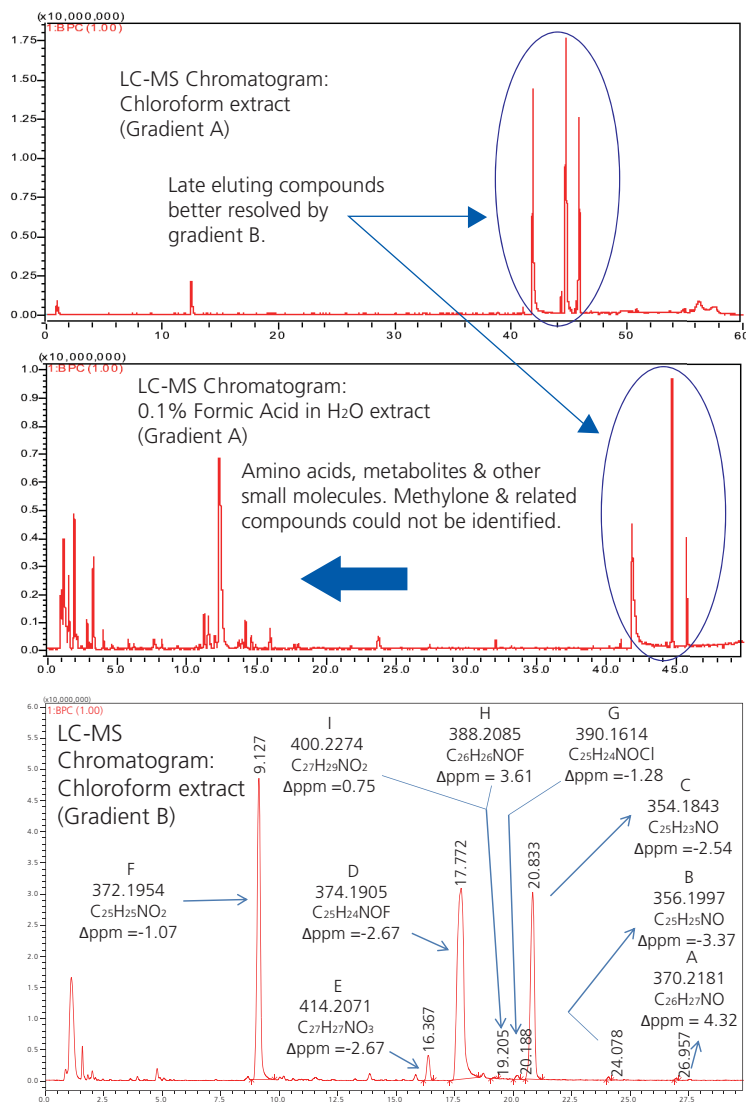


Accurate Mass LC-MS/MS Profiling of Synthetic Cannabinoids

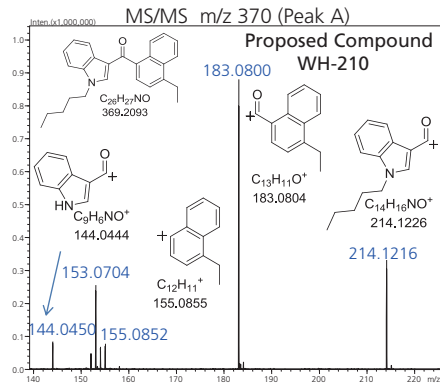

METLIN ID: 64772 NAME: JWH 210 MASS: 369.2093 [View MS/MS] [STRUCTURE]

Frag. m/z	Appm	Intensity	CE	Predicted Ion Type	Predicted Fragment Structure
183.0800	0	100.0	10, 20, 40	[M+H] ⁺	
214.1220	1	33.1	10, 20, 40	[M+H] ⁺	
155.0850	1	73.1	20, 40	[M] ⁺	
144.0440	6	38.4	20, 40	[M+H] ⁺	
153.0900	9	35.3	20, 40	-	No Structure Information available

METLIN: Fragment Ion Metabolite Search
LC Peak A: Proposed Match JWH-210

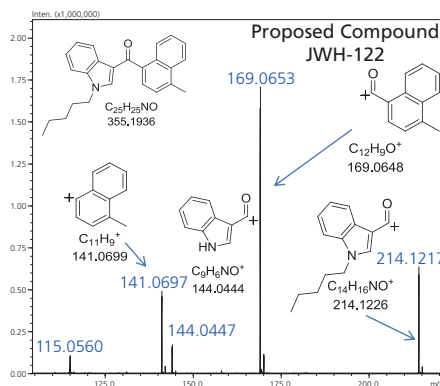


METLIN ID: 64771 NAME: JWH 122 MASS: 355.1936 [View MS/MS] [STRUCTURE]

Frag. m/z	Appm	Intensity	CE	Predicted Ion Type	Predicted Fragment Structure
141.0690	4	100.0	10, 20, 40	[M] ⁺	
214.1200	7	29.1	10, 20, 40	[M+H] ⁺	
169.0630	13	100.0	10, 20, 40	[M+2H] ⁺	
115.0540	17	14.0	20, 40	[M+H] ⁺	

METLIN: Fragment Ion Metabolite Search
LC Peak B: Proposed Match JWH-122

LC Peaks A & B: The proposed structures for LC peaks A and B are known ethyl and methyl naphthyl analogues of synthetic cannabinoid JWH-018. The structures for these two compounds were proposed based on the ability to match accurate mass and MS/MS data with an existing


entry in the METLIN database. The ability to match structures A and B with existing database entries greatly expedited the process of identifying other similar compounds present in the herbal sample extract.

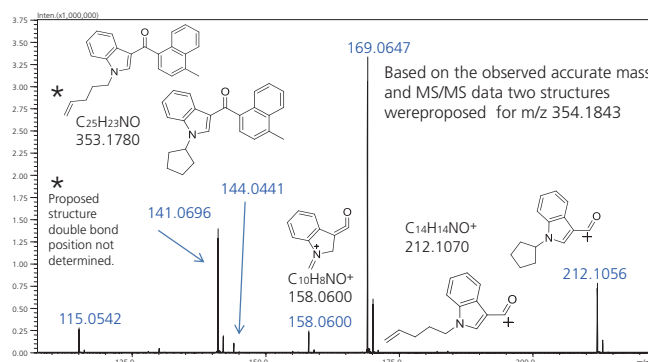

LC Peak C: The MS/MS data for LC Peak C suggested that the compound was a methyl naphthyl compound similar to LC peak B (JWH-122). The MS/MS data further suggested that the difference between compounds B and C was the addition of a double bond or ring ($\Delta = 2$ Da) to the

nitrogen R group. A ChemSpider similarity search of the two structures yielded one compound with a chemical formula of C₂₅H₂₃NO (JWH-150). MS/MS data did not match this structure, however a third structure could be proposed based on this.

Accurate Mass LC-MS/MS Profiling of Synthetic Cannabinoids

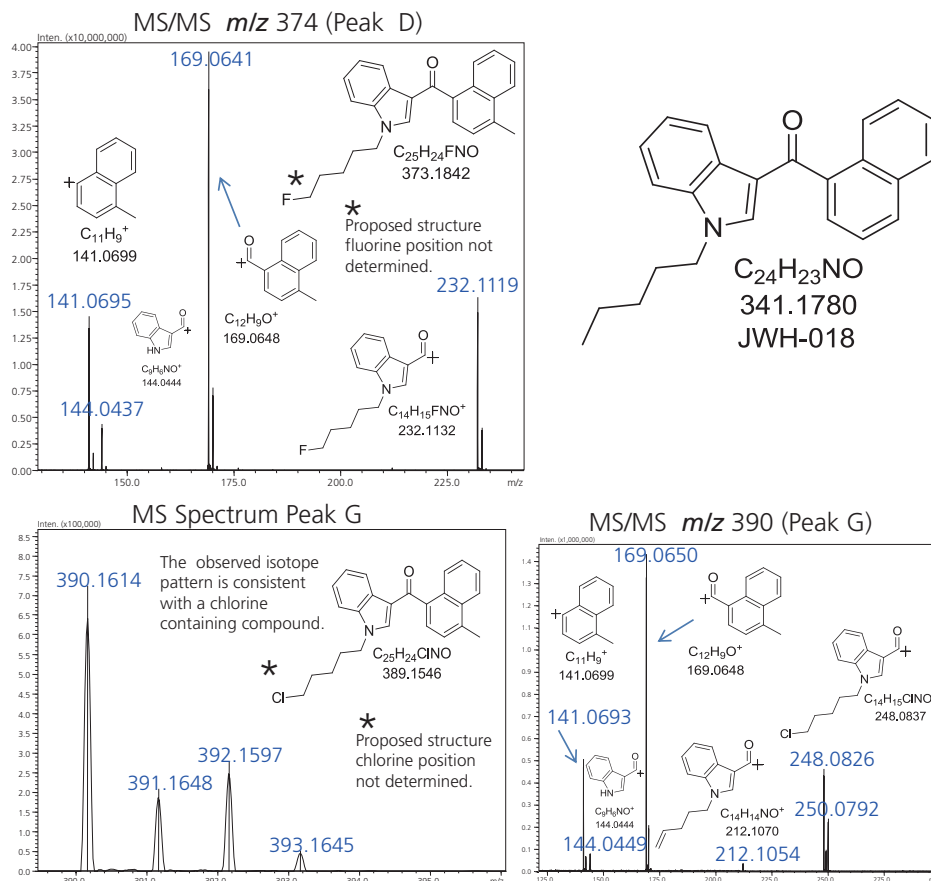

LC Peaks D & G: The proposed structures for compounds D & G are structurally similar to known halogenated analogues of synthetic cannabinoid JWH-018.

Conclusion

From the herbal sample chloroform extract 13 observed compounds had accurate mass and MS/MS data consistent with synthetic cannabinoids.

All compounds gave either a m/z 169 or m/z 183 fragment ion, which suggests that they are either methyl or ethyl naphthyl analogues to JWH-018. Proposed chemical formulas also suggest that some compounds are halogenated.

Of the compounds observed two (JWH-210 & JWH-122) were identified using the METLIN metabolite data base accurate mass & MS/MS ion search. These two compounds were "certified" by the vendor to not be present in this product. It should be noted that this study was not quantitative in nature, and information regarding vendor method detection limits were not available. JWH-210 & JWH 122 may be present at levels less than the PDL of the vendor's method.

From a detection stand point; the results indicate that LC-MS/MS precursor ion scanning on a QqQ or QqToF may be a viable option to screen for classes of synthetic cannabinoids.

Acknowledgements

The authors would like to thank and acknowledge Rice University and the Rice Shared Equipment Authority (SEA) for their support.



Shimadzu Corporation
www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedures.
The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu.
The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.