

FlexiFocus™ slides

Chemically modified sample slides to enhance your MALDI-TOF productivity



Focus your samples in the center of well each time

Perfect for automated workflows

High capacity 96 well slide option for high throughput

Disposable slides to eliminate carryover errors

Enhanced sensitivity performance

FlexiFocus™

... when productivity matters

Product details:

- Single-use metal targets for reproducible centred sample spots
- The slides are composed of a metal substrate which has been chemically modified to 'focus' your samples right where you need them
- The slide consists of an insert with an easy-to-use custom adaptor
- Clearly outlined calibration spots for easy identification
- Flexible option with either a 48 or 96 well slide insert
- Designed for the Shimadzu benchtop MALDI-TOF instruments, AXIMA series*, MALDI-7090* and MALDImini-1**
(*when used with the FlexiMass Precision adaptor (one FlexiFocus adaptor required per slide)
(**MALDImini-1 supports only 48 well slide)

There are 2 separate kits available in the range:

- FlexiFocus-48 well Starter Kit containing adapter and 20 inserts (P/N 224-38626-58)
- FlexiFocus-96 well Starter Kit containing adapter and 20 inserts (P/N 224-38623-58)

Refill items for the kits:

- FlexiFocus-48 well insert, pack of 10 (P/N 224-39906-58)
- FlexiFocus-96 well insert, pack of 10 (P/N 224-39907-58)
- FlexiFocus Adaptor, x 1 (P/N 224-39909-58)

N.B. May require handheld barcode scanner (P/N 224-41045-58)

Kratos Analytical Ltd. is registered to ISO 9001 and ISO 13485.

FlexiFocus, MALDI Solutions, FlexiMass, MALDI-8020 and MALDI-8030 are trademarks of Kratos Analytical Ltd.

Kratos Analytical Ltd. A Shimadzu Group Company


Wharfside, Trafford Wharf Road, Manchester M17 1GP

UK Phone: +44 161 888 4400 Fax. +44 161 888 4402

URL <http://www.shimadzu.com/an/>

Benefits:

- Facilitates automated analysis by producing consistently centred samples spots across the slide
- Improves productivity since analysis performed on center of well each time
- No cleaning prior to use since disposable slide:
 - eliminates carry-over issues between plate cleaning cycles
- Suitable for most types of MALDI samples and matrices
- Provides sample tracking with the 2D barcode label (requires dedicated MALDI Solutions software)
- Higher capacity solution with the FlexiFocus-96 insert
- Enhanced sensitivity performance (*when compared to regular slide under specified experimental conditions)

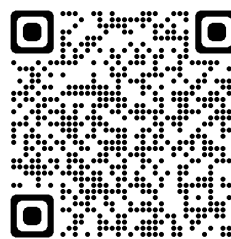

Discover more

For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

Company names, products/service names and logos used in this publication are trademarks and trade names of Shimadzu Corporation, its subsidiaries or its affiliates, whether or not they are used with trademark symbol "TM" or "®".

Third-party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®".

Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

The contents of this publication are provided to you "as is" without warranty of any kind, and are subject to change without notice.

Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication.