


# GC CAPILLARY COLUMN APPLICATIONS: ENVIRONMENTAL

## EPA Method 507 Pesticides/ Herbicides

COLUMN: 007-17, (50% Phenyl)  
Methylpolysiloxane  
25M. x 0.25mm I.D. x  
0.25µm film  
Cat. No.: 007-17-25-0.25F  
Temperature: 80° (6°/min.) - 260°C  
Injector: 230°C  
Detector: MSD  
Carrier Gas: 28 cm/sec., Helium

1. mevinphos (isomer)
2. mevinphos (isomer)
3. demeton
4. phorate
5. disulfoton
6. diazinon and atrazine
7. simazine
8. alachlor
9. fenamophos and bentazone
10. bromocil
11. methyl parathion
12. ametryn
13. malathion
14. ethyl parathion
15. hexazinone
16. ethion
17. azinphos methyl

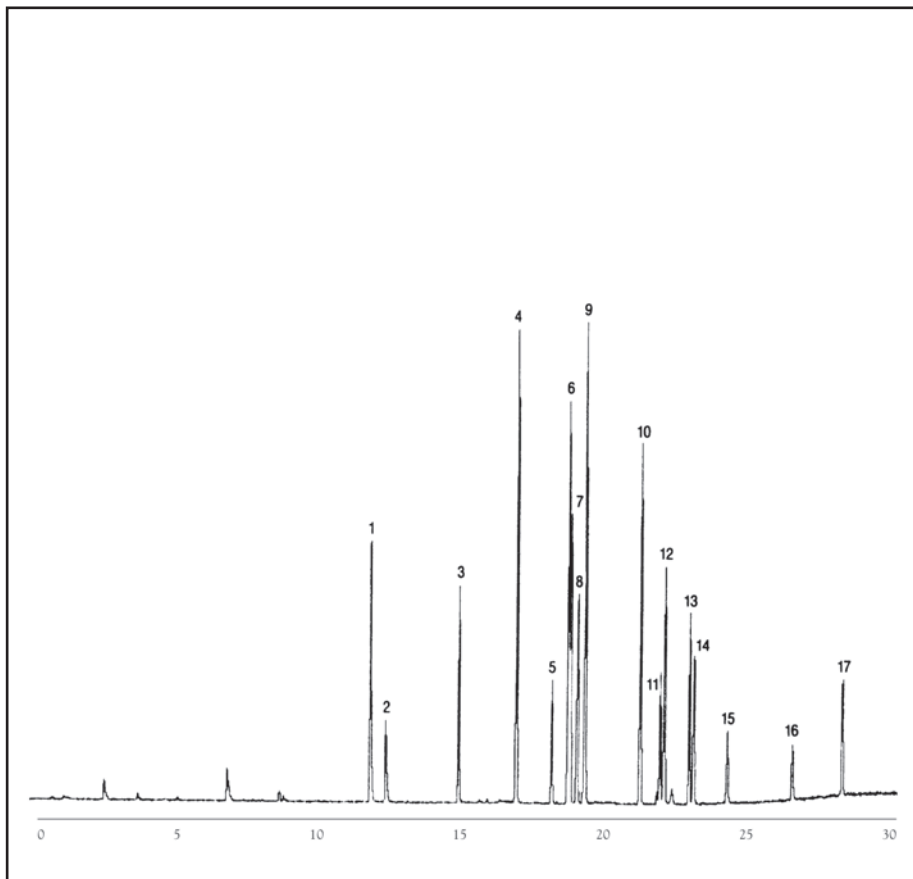

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com)