


GC CAPILLARY COLUMN APPLICATIONS: FAMES

Fame Mix

COLUMN: 007-23,
Bonded Cyanopropyl
Polysiloxane
25M. x 0.25mm I.D. x
0.25µm film
Cat. No.: 007-23-25-0.25F
Temperature: 100° (8°/min.) - 240°C
Injector: 230°C
Detector: 300°C, FID
Carrier Gas: 28 cm/sec., Helium

1. C-14:0
2. C-14:1
3. C-16:0
4. C-16:1 trans
5. C-16:1 cis
6. C-18:0
7. C-18:1 trans
8. C-18:1 cis
9. C-18:2
10. C-20:0
11. C-18:3
12. C-20:1
13. C-20:2
14. C-20:3

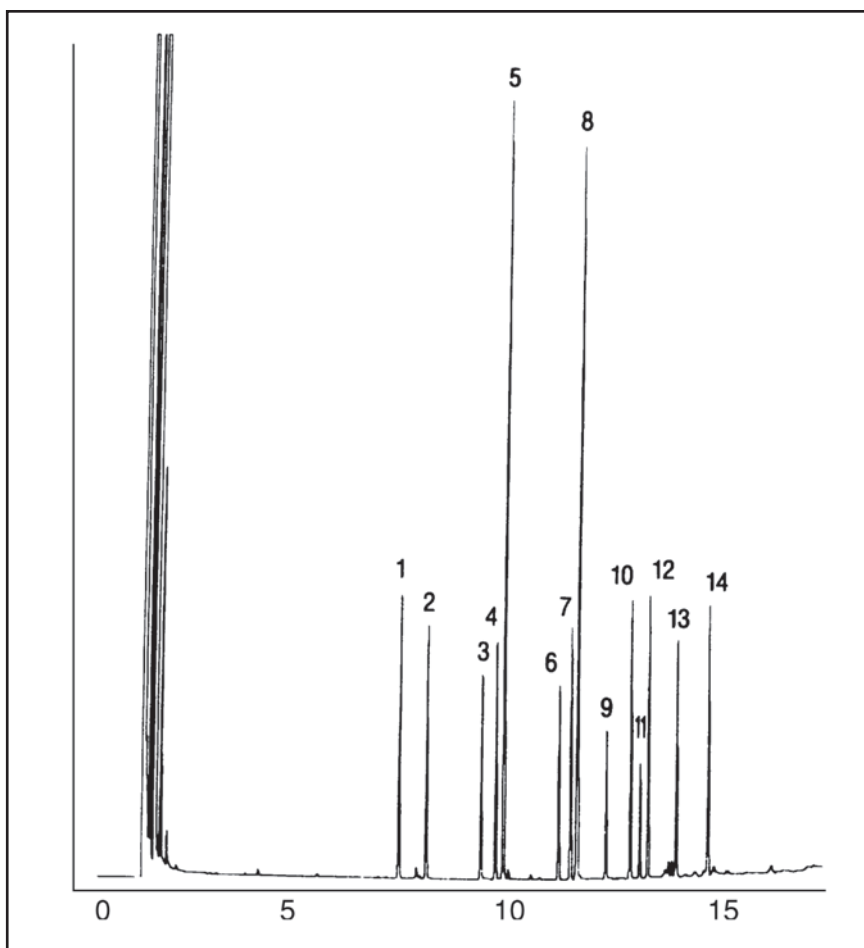

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com