


GC CAPILLARY COLUMN APPLICATIONS: FAMES

Volatile Free Fatty Acids

COLUMN: 007-FFAP,
Bonded FFAP
15M. x 0.32mm I.D. x
0.25µm film
Cat. No.: 007-FFAP-15W-0.25F
Temperature: 110° (10°/min.) - 250°C
Injector: 220°C
Detector: 250°C, FID
Carrier Gas: 26 cm/sec., Helium

1. acetic
2. propionic
3. butyric
4. isovaleric
5. valeric
6. isocaproic
7. caproic
8. enanthic
9. caprylic
10. capric
11. lauric
12. myristic
13. palmitic
14. stearic
15. arachidic

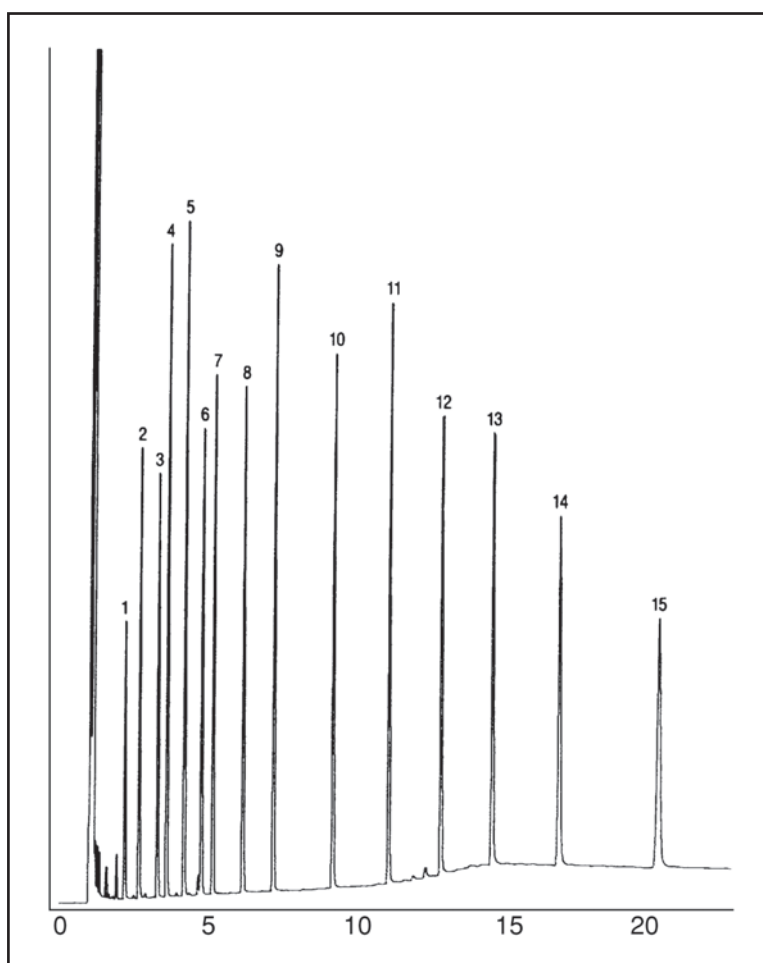

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com